# ASEPTIQUIK® C SERIES CONNECTORS

AseptiQuik® C Connectors provide quick and easy sterile connections, even in non-sterile environments. AseptiQuik® C connectors' "CLICK-PULL-TWIST" design enables users to transfer media easily with less risk of operator error. The connectors' robust design provides reliable performance without the need for clamps, fixtures or tube welders. Biopharmaceutical manufacturers can make sterile connections with the quality and market availability they expect from the leader in single-use connection.



# **SPECIFICATIONS**

**PRESSURE** 

Up to 60 psi, 4.1 bar

**OPERATING TEMPERATURE** 39°F to 104°F (4°C to 40°C)

#### **STERILIZATION**

Standard (blue) version:

**Gamma:** Up to 50 kGy irradiation

High Temperature (white) version:

Gamma: Up to 50 kGy irradiation Autoclave: Up to 266°F (130°C),

60 minutes

#### **TERMINATIONS**

3/8 (9.5mm)" and 1/2" (12.7mm) ID hose barb, and 3/4" sanitary 1-1/2" sanitary

### **MATERIALS**

# **Main Components:**

Polycarbonate (white)

#### Lock Ring:

Polycarbonate (blue)

#### **Pull Tabs**

Polycarbonate (blue, standard version) Polycarbonate (white, HT version)

#### Caps:

Polypropylene (clear)

#### Seals:

Silicone (clear), platinum-cured

## Membrane:

Polyethylene (standard version) Hydrophobic polyethersulfone (HT version),

PTFE strip sticker

WARNING: Pressure, temperature, chemicals, and operating environment can affect the performance of couplings. It is the customer's responsibility to test the suitability of CPC's products in their own application conditions.

#### **FEATURES**

Robust construction —

Integrated lock ring

Membrane pull tabs

CPC Click -

#### **BENEFITS**

Intuitive three-step connection process reduces risk of operator error

Ensure simultaneous and secure removal of both membranes

Repeatable and reliable performance with no additional hardware required

Secures final connection preventing disassembly

Audible confirmation of completed assembly steps

#### TYPICAL FLOW RATE:

Cv Value Range: 4-14.4



Cv values represent the approximate expected flow rate in gallons per minute of water at room temperature for a 1 PSI pressure drop. The flow is generally constrained by the smallest diameter, which in some cases will be the termination diameter and not the Nominal Flow Path.

Scan code to visit webpage



cpcworldwide.com/AseptiQuik-C

# **DID YOU KNOW**

Did you know that the AseptiQuik G Series Connectors also have 3/8", 1/2", and 3/4" sanitary connections?



COLDER PRODUCTS COMPANY

D.S.A.

+1 (651) 645-0091

X: +1 (651) 645-5404 DLL FREE: +1 (800) 519-7633 MAIL: info@cpcworldwide.com COLDER PRODUCTS COMPANY GMBH
Germany

cocgmbh@cocworldwide.com

DOVER (SHANGHAI) INDUSTRIAL CO., LTD ShangHai. China

PHONE: +86 21 2411 2666 TOLL FREE: +86 400 990 1978

E-MAIL: asiapacific@cpcworldwide.co

# **ASEPTIQUIK® C SERIES DIMENSIONS**

# **COUPLING BODY** Polycarbonate - with Blue Pull Tabs



TERMINATION	METRIC EQ.	PART NO.	A	В	<b>©</b>
3/8" HOSE BARB	9.5 mm	AQC17006	2.36" (59.9mm)	2.74" (69.6mm)	0.65" (16.5mm)
1/2" HOSE BARB	12.7 mm	AQC17008	2.36" (59.9mm)	2.95" (74.9mm)	0.89" (22.6mm)
3/4" SANITARY	19.1 mm	AQC33012	2.36" (59.9mm)	2.73" (69.3mm)	0.70" (17.8mm)

# **COUPLING BODY** Polycarbonate - *HT with White Pull Tabs*



TERMINATION	METRIC EQ.	PART NO.	A	В	0
3/8" HOSE BARB	9.5 mm	AQC17006HT	2.36" (59.9mm)	2.74" (69.6mm)	0.65" (16.5mm)
1/2" HOSE BARB	12.7 mm	AQC17008HT	2.36" (59.9mm)	2.95" (74.9mm)	0.89" (22.6mm)
3/4" SANITARY	19.1 mm	AQC33012HT	2.36" (59.9mm)	2.73" (69.3mm)	0.70" (17.8mm)

# **COUPLING INSERT** Polycarbonate - with Blue Pull Tabs




TERMINATION	METRIC EQ.	PART NO.	A	В	C
3/8" HOSE BARB	9.5 mm	AQC22006	2.50" (63.5mm)	2.90" (73.7mm)	0.65" (16.5mm)
1/2" HOSE BARB	12.7 mm	AQC22008	2.50" (63.5mm)	2.99" (75.9mm)	0.89" (22.6mm)
3/4" SANITARY	19.1 mm	AQC44012	2.50" (63.5mm)	2.80" (71.1mm)	0.70" (17.8mm)


# **COUPLING INSERT** Polycarbonate - HT with White Pull Tabs



TERMINATION	METRIC EQ.	PART NO.	A	<b>B</b>	G
3/8" HOSE BARB	9.5 mm	AQC22006HT	2.50" (63.5mm)	2.90" (73.7mm)	0.65" (16.5mm)
1/2" HOSE BARB	12.7 mm	AQC22008HT	2.50" (63.5mm)	2.99" (75.9mm)	0.89" (22.6mm)
3/4" SANITARY	19.1 mm	AQC44012HT	2.50" (63.5mm)	2.80" (71.1mm)	0.70" (17.8mm)

# PRODUCT DIMENSIONS

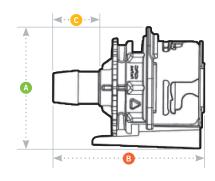

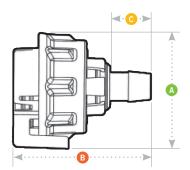


A = Height/Diameter

B = Total Length

c = Hose Barb Length

# **AQC ASSEMBLY PROCEDURE**

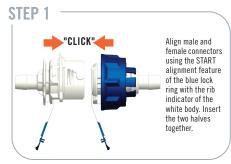



Align male and female connectors using the START alignment feature of the blue lock ring with the rib indicator of the white body. Insert the two halves together.



Align male and female connectors using the START alignment feature of the blue lock ring with the rib indicator of the white body. Insert the two halves together.

VARRANTY: All sales are subject to Colder Products Company's limited express warranty set forth in the CPC catalog. Contact your local distributor or CPC Customer Service for warranty provision.

Warring: Due to the wide variety of possible fluid media and operating conditions, unintended consequences may result from the use of this product, all of which are beyond the control of CPC. It is the user's responsibility to carefully determine and test for compatibility for use with their application. All such risks shall be assumed by the buyer.

COPYRIGHT © 2022 BY COLDER PRODUCTS COMPANY