

- Made From USP Class VI Material
- Assembled in Cleanroom 7 Cleanroom
- Manufactured according to ISO 13485 standards
- All Connectors Individually Double Bagged with Labels on the inner bag
- Available Sterile or Non-Sterile
- Available with 3/4" or 1-1/2" Sanitary Termination of Steam Thru Connector and AseptiQuik G or AseptiQuik S versions of Sterile Connectors

• **Manufactured in USA**




**Foxx Life Sciences New Steam-Thru AseptiQuik Series Connector** allows a quick and easy sterile connection between biopharmaceutical processing equipment and single-use systems. **AseptiQuik® STC** Connectors integrate the **AseptiQuik** sterile connector and the **Steam-Thru® II SIP** connector, giving manufacturers flexibility between hybrid stainless steel and single-use processing equipment. The **AseptiQuik STC** connector combines either an **AseptiQuik G** or **AseptiQuik S** connector with a **Steam-Thru II** connection that can be mounted directly to stainless steel processing equipment via a sanitary termination.

The Union of the two connectors into a single unit is through a Foxx sanitary Tri-clamp allows **AseptiQuik** sterile connection to be steamed onto the stainless-steel equipment via SIP. After the SIP cycle, a wide range of Single-Use Systems can be connected with **AseptiQuik G** or **AseptiQuik S** Sterile connector. The SIP process can be done in advance, allowing a quick and easy sterile connection to the **AseptiQuik** half without waiting 30-60 minutes for SIP before harvest. The **AseptiQuik STC** Connector connects stainless steel processing equipment via 3/4" or 1-1/2" Sanitary Termination.

- 1 EZBio® Clamp
  - Can withstand 10 bars pressure
- 2 AseptiQuik® Available in:
  - Standard and HT version
- 3 Steam-Thru® Available in:
  - 3/4" x 3/4" Sanitary
  - 3/4" x 1 1/2" Sanitary



ISO 13485

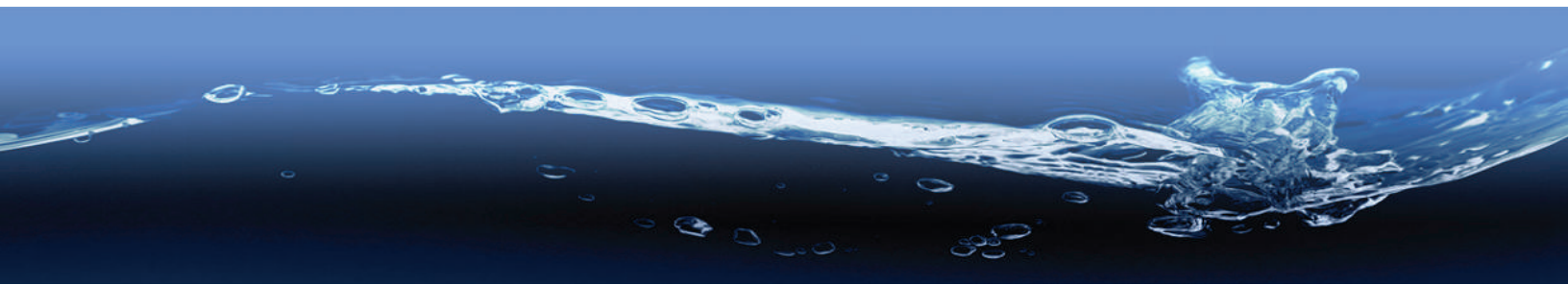

	Part Number	Description	AseptiQuik® Body	Sanitary Termination	Autoclavable	Qty/Pack
Sterile	446-301C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik G Series Connector	G-Series	3/4" x 1-1/2"	No	5/CS
	446-301D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik G Series Connector	G-Series	3/4" x 3/4"	No	5/CS
	446-341C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik S Series Connector	S-Series	3/4" x 1-1/2"	No	5/CS
	446-341D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik S Series Connector	S-Series	3/4" x 3/4"	No	5/CS
	446-311C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT G Series Connector	HT G-Series	3/4" x 1-1/2"	Yes	5/CS
	446-311D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT G Series Connector	HT G-Series	3/4" x 3/4"	Yes	5/CS
	446-351C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT S Series Connector	HT S-Series	3/4" x 1-1/2"	Yes	5/CS
	446-351D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT S Series Connector	HT S-Series	3/4" x 3/4"	Yes	5/CS
Non-Sterile	446-301C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik G Series Connector	G-Series	3/4" x 1-1/2"	No	5/CS
	446-301D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik G Series Connector	G-Series	3/4" x 3/4"	No	5/CS
	446-341C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik S Series Connector	S-Series	3/4" x 1-1/2"	No	5/CS
	446-341D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik S Series Connector	S-Series	3/4" x 3/4"	No	5/CS
	446-311C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT G Series Connector	HT G-Series	3/4" x 1-1/2"	Yes	5/CS
	446-311D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT G Series Connector	HT G-Series	3/4" x 3/4"	Yes	5/CS
	446-351C-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT S Series Connector	HT S-Series	3/4" x 1-1/2"	Yes	5/CS
	446-351D-FLS	EZBio® Custom MTO SU Connector Assembly Steam-Thru w/ Aseptiquik HT S Series Connector	HT S-Series	3/4" x 3/4"	Yes	5/CS