


EZwaste® HPLC Caps

www.foxxlifesciences.com



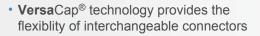

- To prevent spills and vapors from leaking, tube fittings secure tubes to container to create a close-loop system.
- Plugs included for all unused ports

C MACHINES

• Depending on application, carbon exhaust filter life cycle may vary

Solvent V

Foxx Life Sciences' **EZ**Waste® HPLC Caps are designed to be an affordable solution for liquid waste disposal. Each system is manufactured to combine with any 38-430, GL45, 51S, 53, S70 or 83B threaded container and safely collect waste from one or more HPLC systems. These assemblies feature an activated carbon exhaust filter which traps volatile organic compounds form being released into the air. This prevents chemicals from being inhaled by lab workers and provides an easy mechanism of disposal for vapor wastes. The **EZ**Waste® HPLC Caps feature a closed loop design which protects against leaks and spills, improving overall lab safety.

Each cap assembly includes Delrin or polypropylene fittings with corresponding plugs for unused ports. Furthermore, each assembly will come with a **Versa**Cap® 38-430, GL45, 51S, 53, S70 or 83B with a corresponding manifold and a carbon exhaust filter.

- 1 Carbon exhaust filter traps* solvent vapor
- 2 Polypropylene 1/2" OD compression tube fitting
- 3 1/8" and/or 3/8"Polypropelyne Hose Barb Fittings
- 4 Delrin 1/8" and/or 1/16" OD tube fittings with corresponding Ferrule(s)
- 5 EZWaste® Universal Stackable Manifold with 6 Ports for Either 1/16" or 1/8" OD Tubing, 3 Ports for Either 1/8" or 3/8" Hose Barbs, 3 Ports For 1/2" OD Tubing and Matching Plugs for Each Port
- 6 VersaCap® 38-430, Open Top for Use with Machines Manifolds, White PP

Made in the USA

ISO 13485:2016



EZwaste® HPLC Caps

www.foxxlifesciences.com



	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
38-430	EZwaste® Cap Assembly, 3A Ports	Flat	Yes	330-0101-FLS
	EZwaste® Cap Assembly, 6A Ports, 3C Ports, 3D Ports	Universal	Yes	330-0103-FLS
	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
GL45	EZwaste® Cap Assembly, 6A Ports	Flat	Yes	330-0D11-FLS
	EZwaste® Cap Assembly, 4A Ports, 2B Ports	Flat	Yes	330-0D21-FLS
	EZwaste® Cap Assembly,6A Ports, 3C Ports, 3D Ports	Universal	Yes	330-0D31-FLS
	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
51	EZwaste® Cap Assembly, 6A Ports	Flat	Yes	330-0A53-FLS
	EZwaste® Cap Assembly, 4A Ports, 2B Ports	Flat	Yes	330-0A51-FLS
	EZwaste® Cap Assembly, 6A Ports, 3C Ports, 3D Ports	Universal	Yes	330-0A52-FLS
	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
53	EZwaste® Cap Assembly, 6A Ports	Flat	Yes	330-0212-FLS
	EZwaste® Cap Assembly, 4A Ports, 2B Ports	Flat	Yes	330-0213-FLS
	EZwaste® Cap Assembly, 6A Ports, 3C Ports, 3D Ports	Universal	Yes	330-0214-FLS
	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
S70	EZwaste® Cap Assembly, 6A Ports, 3C Ports	Universal	Yes	330-0B60-FLS
	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
83B	EZwaste® Cap Assembly, 6A Ports	Flat	Yes	330-0441-FLS
	EZwaste® Cap Assembly, 4A Ports, 2B Ports	Flat	Yes	330-0457-FLS
	EZwaste® Cap Assembly, 6A Ports, 3C Ports, 3D Ports	Universal	Yes	330-0458-FLS
	DESCRIPTION	MANIFOLD TYPE	EXHAUST FILTER	PART NUMBER
N/A	EZwaste® Manifold, 6A Ports, 3C Ports, 3D Ports	Universal	Yes	330-0E32-FLS