


PharmaFluor® PFA Compression Fittings

Saint-Gobain manufactures PFA fittings for handling corrosive or ultra-pure fluids under contamination-free conditions. The PharmaFluor® compression fitting line offers PFA purity combined with the robust performance of compression style fittings. These fittings are ideal for gas or pneumatic connections in applications that do not require flared connections. Once tubing has been cut to the desired length, PharmaFluor® fittings require no tools for installation. The fittings are available in 1/8", 1/4", 3/8", 1/2", 3/4" and 1" sizes.

We use 100% virgin HP PFA to give you the ultimate in purity and the minimum in contamination. Each PharmaFluor® PFA fitting is carefully inspected to ensure quality and reliability in your application and throughout the life of the product. Fluid handling applications in the biomedical, chemical, petrochemical, food processing and semiconductor industries benefit greatly from the use of these fittings.

Whether it be a simple fluid handling system or a complex laboratory layout, PharmaFluor® PFA fittings save precious time and money, and provide all the reliability needed to facilitate successful outcomes.

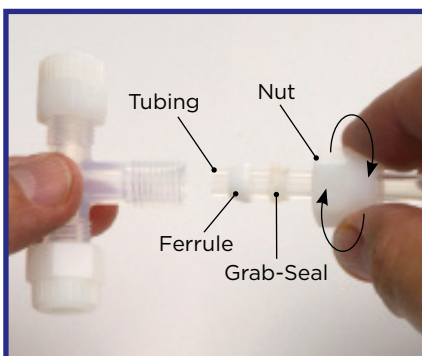

Fitting Installation

1. Cut tubing to desired length.
2. Loosen nut, insert tube.
(your pre-assembled fitting should be ready for use)
3. Tighten nut finger-tight.
4. Your fitting is now ready for use.

No tools are required, but if you decide to use a mechanical device to tighten it, rotate the nut no more than one turn past finger-tight.


Features / Benefits

- Simple design for ease of installation
- No special tools required
- Accommodates temperatures up to 212°F (100°C)

Materials of Construction


- Body: high purity PFA
- Gripper: PVDF
- Ferrule: PTFE
- Nut: ETFE

Union Adapter Male
Male NPT Threads

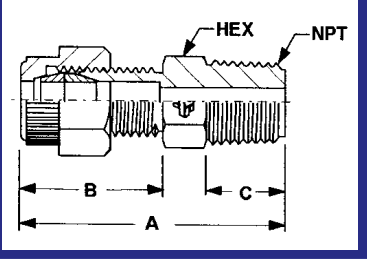


Part Number	O.D Tube Size Inches (mm)	Male NPT Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	Hex Flat Inches (mm)
G1077067	1/8 (3.17)	1/8 (3.17)	0.09 (2.29)	1.40 (35.56)	0.73 (18.54)	0.40 (10.16)	0.45 (11.43)
G1077054	1/8 (3.17)	1/4 (6.35)	0.09 (2.29)	1.50 (38.10)	0.73 (18.54)	0.50 (12.71)	0.50 (12.70)
G1077053	1/4 (6.35)	1/8 (3.17)	0.19 (4.83)	1.47 (37.34)	0.80 (20.32)	0.40 (10.16)	0.44 (11.18)
G1077002	1/4 (6.35)	1/4 (6.35)	0.19 (4.83)	1.58 (40.13)	0.81 (20.57)	0.50 (12.70)	0.50 (12.70)
G1077080	1/4 (6.35)	3/8 (9.52)	0.19 (4.83)	1.67 (42.42)	0.80 (20.32)	0.50 (12.70)	0.68 (17.27)
G1077081	1/4 (6.35)	1/2 (12.70)	0.19 (4.83)	1.80 (45.72)	0.80 (20.32)	0.63 (16.00)	0.87 (22.10)
G1077382	3/8 (9.52)	1/8 (3.17)	0.19 (4.83)	1.61 (40.89)	0.84 (21.34)	0.40 (10.16)	0.68 (17.27)
G1077003	3/8 (9.52)	1/4 (6.35)	0.31 (7.87)	1.71 (43.43)	0.84 (21.34)	0.50 (12.71)	0.68 (17.27)
G1077004	3/8 (9.52)	3/8 (9.52)	0.31 (7.87)	1.71 (43.43)	0.84 (21.34)	0.50 (12.71)	0.68 (17.27)
G1077051	3/8 (9.52)	1/2 (12.70)	0.31 (7.87)	1.84 (46.74)	0.84 (21.34)	0.63 (16.00)	0.87 (22.10)
G1077050	1/2 (12.70)	1/4 (6.35)	0.29 (7.37)	1.81 (45.97)	0.94 (23.88)	0.50 (12.71)	0.87 (22.10)
G1077005	1/2 (12.70)	3/8 (9.52)	0.42 (10.67)	1.81 (45.97)	0.94 (23.88)	0.50 (12.71)	0.87 (22.10)
G1077001	1/2 (12.70)	1/2 (12.70)	0.46 (11.68)	1.95 (49.53)	0.94 (23.88)	0.64 (16.26)	0.87 (22.10)
G1077104	1/2 (12.70)	3/4 (19.05)	0.46 (11.68)	2.18 (55.37)	0.88 (22.35)	0.80 (20.32)	1.12 (28.45)

Union Adapter Female
Female NPT Threads

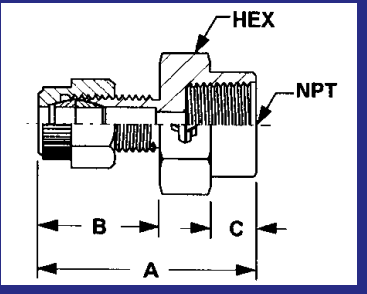


Part Number	O.D Tube Size Inches (mm)	Female NPT Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	Hex Flat Inches (mm)
G1077084	1/4 (6.35)	1/8 (3.17)	0.19 (4.83)	1.38 (35.05)	0.80 (20.32)	0.31 (7.87)	0.56 (14.22)
G1077055	1/4 (6.35)	1/4 (6.35)	0.19 (4.83)	1.52 (38.61)	0.80 (20.32)	0.34 (8.64)	0.94 (23.88)
G1077085	1/4 (6.35)	3/8 (9.52)	0.19 (4.83)	1.48 (37.59)	0.80 (20.32)	0.31 (7.87)	1.00 (25.40)
G1077086	1/4 (6.35)	1/2 (12.70)	0.19 (4.83)	1.57 (39.88)	0.80 (20.32)	0.40 (10.16)	1.13 (28.70)
G1077377	3/8 (9.52)	1/8 (3.17)	0.25 (6.35)	1.52 (38.61)	0.84 (21.34)	0.31 (7.87)	0.69 (17.53)
G1077056	3/8 (9.52)	1/4 (6.35)	0.31 (7.87)	1.53 (38.86)	0.84 (21.34)	0.32 (8.13)	0.94 (23.88)
G1077057	3/8 (9.52)	3/8 (9.52)	0.31 (7.87)	1.52 (38.61)	0.84 (21.34)	0.31 (7.87)	1.00 (25.40)
G1077087	3/8 (9.52)	1/2 (12.70)	0.31 (7.87)	1.61 (40.89)	0.84 (21.34)	0.40 (10.16)	1.13 (28.70)
G1077088	1/2 (12.70)	1/4 (6.35)	0.31 (7.87)	1.61 (40.89)	0.94 (23.88)	0.30 (7.62)	0.88 (22.35)
G1077089	1/2 (12.70)	3/8 (9.52)	0.45 (11.43)	1.63 (41.40)	0.94 (23.88)	0.32 (8.13)	1.13 (28.70)
G1077090	1/2 (12.70)	1/2 (12.70)	0.46 (11.68)	1.73 (43.94)	0.94 (23.88)	0.42 (10.67)	1.13 (28.70)
G1077373	1/2 (12.70)	3/4 (19.05)	0.46 (11.68)	1.82 (46.23)	0.94 (23.88)	0.38 (9.63)	1.38 (35.05)

Elbow Adapter Male
Male NPT Threads

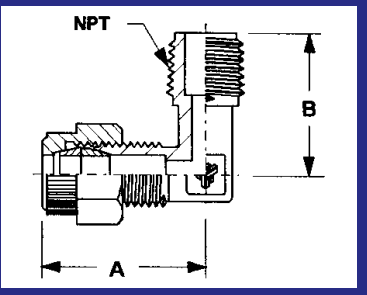


Part Number	O.D Tube Size Inches (mm)	Male NPT Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)
G1077072	1/8 (3.17)	1/8 (3.17)	0.09 (2.29)	0.73 (18.54)	0.72 (18.29)
G1077091	1/8 (3.17)	1/4 (6.35)	0.09 (2.29)	0.99 (25.15)	0.78 (19.81)
G1077059	1/4 (6.35)	1/8 (3.17)	0.19 (4.83)	1.05 (26.67)	0.71 (18.03)
G1077007	1/4 (6.35)	1/4 (6.35)	0.19 (4.83)	1.05 (26.67)	0.97 (24.64)
G1077321	1/4 (6.35)	3/8 (9.52)	0.19 (4.83)	1.16 (29.46)	0.85 (21.59)
G1077323	3/8 (9.52)	1/8 (3.17)	0.19 (4.83)	1.19 (30.23)	0.85 (21.59)
G1077008	3/8 (9.52)	1/4 (6.35)	0.31 (7.87)	1.19 (30.23)	0.95 (24.13)
G1077009	3/8 (9.52)	3/8 (9.52)	0.31 (7.87)	1.18 (29.97)	1.05 (26.67)
G1077094	3/8 (9.52)	1/2 (12.70)	0.31 (7.87)	1.28 (32.51)	1.13 (28.70)
G1077058	1/2 (12.70)	1/4 (6.35)	0.31 (7.87)	1.28 (32.51)	1.13 (28.70)
G1077010	1/2 (12.70)	3/8 (9.52)	0.42 (10.67)	1.28 (32.51)	1.13 (28.70)
G1077012	1/2 (12.70)	1/2 (12.70)	0.46 (11.68)	1.36 (34.54)	1.13 (28.70)
G1077105	1/2 (12.70)	3/4 (19.05)	0.46 (11.68)	1.57 (39.88)	1.50 (38.10)

Elbow Adapter Female
Female NPT Threads

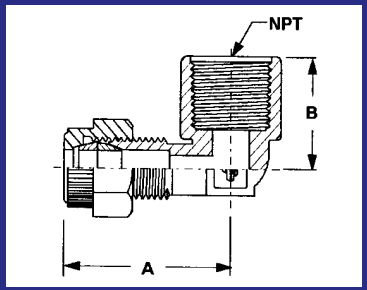


Part Number	O.D Tube Size Inches (mm)	Female NPT Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)
G1077317	1/4 (6.35)	1/4 (6.35)	0.09 (2.29)	1.27 (32.26)	0.80 (20.32)
G1077318	3/8 (9.52)	1/4 (6.35)	0.31 (7.87)	1.32 (33.53)	1.00 (25.40)
G1077095	3/8 (9.52)	3/8 (9.52)	0.31 (7.87)	1.40 (35.56)	1.00 (25.40)
G1077316	3/8 (9.52)	1/2 (12.70)	0.31 (7.87)	1.59 (40.39)	1.12 (28.45)
G1077322	1/2 (12.70)	3/8 (9.52)	0.44 (11.18)	1.50 (38.10)	1.12 (28.45)
G1077096	1/2 (12.70)	1/2 (12.70)	0.40 (10.16)	1.63 (41.40)	1.25 (31.75)

Male Branch Tee
Male NPT Threads

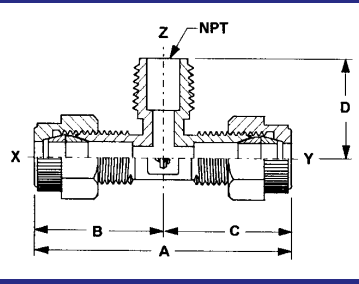


Part Number	O.D Tube Size (X) Inches (mm)	O.D Tube Size (Y) Inches (mm)	Male NPT Size (Z) Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	D Inches (mm)
G1077097	1/8 (3.17)	1/8 (3.17)	1/8 (3.17)	0.09 (2.29)	1.86 (47.24)	0.93 (23.62)	0.93 (23.62)	0.78 (19.81)
G1077047	1/4 (6.35)	1/4 (6.35)	1/4 (6.35)	0.19 (4.83)	2.10 (53.34)	1.05 (26.67)	1.05 (26.67)	0.93 (23.62)
G1077049	3/8 (9.52)	3/8 (9.52)	3/8 (9.52)	0.31 (7.87)	2.38 (60.45)	1.19 (30.23)	1.19 (30.23)	0.94 (23.88)
G1077100	1/2 (12.70)	1/2 (12.70)	1/2 (12.70)	0.44 (11.18)	2.78 (70.61)	1.39 (35.31)	1.39 (35.31)	1.23 (31.24)

Male Run Tee
Male NPT Threads

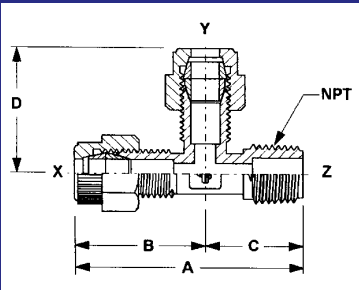


Part Number	O.D Tube Size (X) Inches (mm)	O.D Tube Size (Y) Inches (mm)	Male NPT Size (Z) Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	D Inches (mm)
G1077101	1/8 (3.17)	1/8 (3.17)	1/8 (3.17)	0.09 (2.29)	1.70 (43.18)	0.92 (23.37)	0.78 (19.81)	0.92 (23.37)
G1077039	1/4 (6.35)	1/4 (6.35)	1/4 (6.35)	0.19 (4.83)	1.90 (48.26)	1.05 (26.67)	0.85 (21.59)	1.05 (26.67)
G1077041	3/8 (9.52)	3/8 (9.52)	3/8 (9.52)	0.31 (7.87)	2.13 (54.10)	1.19 (30.23)	0.94 (23.88)	1.19 (30.23)
G1077338	1/2 (12.70)	1/2 (12.70)	1/2 (12.70)	0.44 (11.18)	2.62 (66.55)	1.39 (35.31)	1.23 (31.24)	1.39 (35.31)

Straight Union

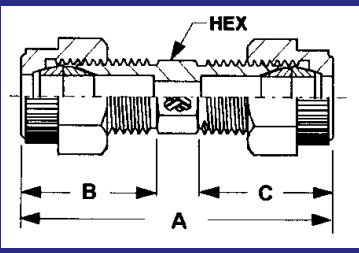


Part Number	O.D Tube Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	Hex Flat Inches (mm)
G1077070	1/8 (3.17)	0.09 (2.29)	1.84 (46.74)	0.73 (18.54)	0.73 (18.54)	0.38 (9.65)
G1077019	1/4 (6.35)	0.19 (4.83)	1.87 (47.50)	0.80 (20.32)	0.80 (20.32)	0.38 (9.65)
G1077020	3/8 (9.52)	0.31 (7.87)	2.05 (52.07)	0.84 (21.34)	0.84 (21.34)	0.68 (17.27)
G1077021	1/2 (12.70)	0.44 (11.18)	2.23 (56.64)	0.93 (23.62)	0.93 (23.62)	0.87 (22.10)

Union Elbow

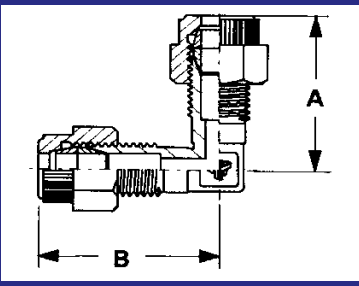


Part Number	O.D Tube Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)
G1077071	1/8 (3.17)	0.09 (2.29)	0.94 (23.88)	0.94 (23.88)
G1077023	1/4 (6.35)	0.19 (4.83)	1.19 (30.23)	1.19 (30.23)
G1077024	3/8 (9.52)	0.31 (7.87)	1.31 (33.27)	1.31 (33.27)
G1077025	1/2 (12.70)	0.44 (11.18)	1.47 (37.34)	1.47 (37.34)

Union Tee

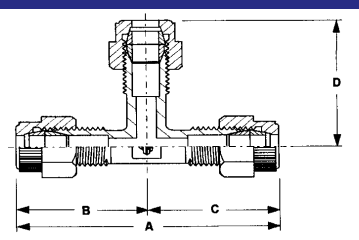


Part Number	O.D Tube Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	D Inches (mm)
G1077066	1/8 (3.17)	0.09 (2.29)	1.86 (47.24)	0.93 (23.62)	0.93 (23.62)	0.93 (23.62)
G1077027	1/4 (6.35)	0.19 (4.83)	2.28 (57.91)	1.14 (28.96)	1.14 (28.96)	1.35 (34.29)
G1077028	3/8 (9.52)	0.31 (7.87)	2.38 (60.45)	1.19 (30.23)	1.19 (30.23)	1.19 (30.23)
G1077029	1/2 (12.70)	0.44 (11.18)	3.16 (80.26)	1.58 (40.13)	1.58 (40.13)	1.58 (40.13)

Straight Union Reducer

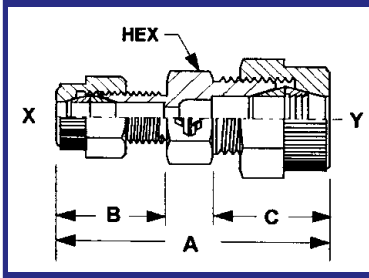


Part Number	O.D Tube Size (X) Inches (mm)	O.D Tube Size (Y) Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	D Hex Flat
G1077060	1/8 (3.17)	1/4 (6.35)	0.09 (2.29)	1.80 (45.72)	0.73 (18.54)	0.81 (20.57)	0.38
G1077061	1/4 (6.35)	3/8 (9.52)	0.19 (4.83)	2.01 (51.05)	0.80 (20.32)	0.84 (21.34)	0.68
G1077062	1/4 (6.35)	1/2 (12.70)	0.19 (4.83)	2.11 (53.59)	0.80 (20.32)	0.94 (23.88)	0.93
G1077063	3/8 (9.52)	1/2 (12.70)	0.31 (7.87)	2.15 (54.61)	0.84 (21.34)	0.94 (23.88)	0.87
G1077103	1/2 (12.70)	3/4 (19.05)	0.70 (17.78)	2.51 (63.75)	0.94 (23.88)	1.07 (27.18)	1.375

Union Elbow Reducer

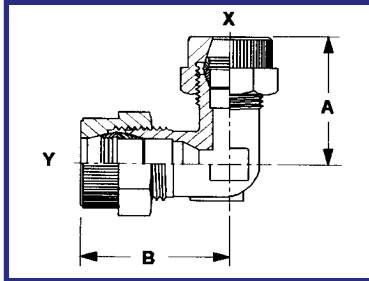


Part Number	O.D Tube Size (X) Inches (mm)	O.D Tube Size (Y) Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)
G1077388	1/8 (3.17)	1/4 (6.35)	0.09 (2.29)	1.00 (25.40)	1.08 (27.43)
G1077387	1/4 (6.35)	3/8 (9.52)	0.19 (4.83)	1.16 (29.46)	1.25 (31.75)
G1077390	3/8 (9.52)	1/2 (12.70)	0.31 (7.87)	1.36 (34.54)	1.44 (36.58)
G1077386	1/4 (6.35)	1/2 (12.70)	0.25 (6.35)	1.30 (33.02)	1.44 (36.58)
G1077384	1/2 (12.70)	3/4 (19.05)	0.44 (11.18)	1.74 (44.20)	1.81 (45.97)

Union Tee Reducer

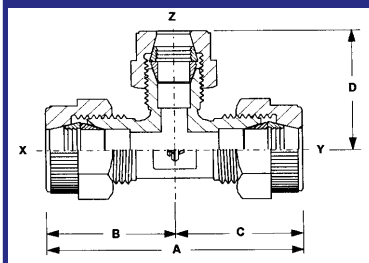

Part Number	O.D. Tube Size (X) Inches (mm)	O.D. Tube Size (Y) Inches (mm)	O.D. Tube Size (Z) Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	D Inches (mm)
G1077395	3/8 (9.52)	3/8 (9.52)	1/4 (6.35)	0.19 (4.83)	2.40 (60.96)	1.20 (30.48)	1.20 (30.48)	1.16 (29.46)
G1077392	1/2 (12.70)	1/2 (12.70)	1/4 (6.35)	0.19 (4.83)	2.80 (71.12)	1.40 (35.56)	1.40 (35.56)	1.32 (33.53)
G1077393	1/2 (12.70)	1/2 (12.70)	3/8 (9.52)	0.31 (7.87)	2.80 (71.12)	1.40 (35.56)	1.40 (35.56)	1.36 (34.54)

Straight Union Panel

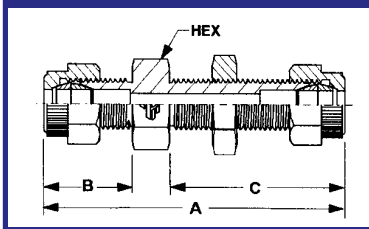

Part Number	O.D Tube Size Inches (mm)	Orifice Inches (mm)	A Inches (mm)	B Inches (mm)	C Inches (mm)	Hex Flat Inches (mm)
G1077068	1/8 (3.17)	0.09 (2.29)	2.80 (71.12)	0.73 (18.54)	1.69 (42.93)	0.58 (14.73)
G1077013	1/4 (6.35)	0.19 (4.83)	2.92 (74.17)	0.80 (20.32)	1.75 (44.45)	0.68 (17.27)
G1077014	3/8 (9.52)	0.31 (7.87)	3.04 (77.22)	0.84 (21.34)	1.83 (46.48)	0.94 (23.88)
G1077015	1/2 (12.70)	0.44 (11.18)	3.16 (80.26)	0.93 (23.62)	1.86 (47.24)	1.12 (28.45)

PharmaFluor® is a registered trademark.

Uncontrolled Document - for the controlled version of this document please visit www.biopharm.saint-gobain.com

Contact us today for:

Consultations • Samples • Quotes • Orders • Technical Service

For registered access to Saint-Gobain Product Validation Summaries, [CLICK HERE](#)

IMPORTANT: It is the user's responsibility to ensure the suitability and safety of Saint-Gobain Life Sciences products for all intended uses and that the materials to be used comply with all applicable medical regulatory requirements. Saint-Gobain Life Sciences assumes no responsibility for any product failures that occur due to misuse of the materials it provides arising out of the design, fabrication or application of the products into which the materials are incorporated.

WARRANTY: For a period of 12 months from the date of first sale, Saint-Gobain Life Sciences warrants this product to be free of defects in materials and workmanship. Our only obligation will be to replace any portion proving defective, or at our option, to refund the purchase price there of.

SAINT-GOBAIN LIFE SCIENCES DISCLAIMS ALL IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.



Saint-Gobain Life Sciences

210 Harmony Road
Mickelton, NJ 08056
Tel: (856) 423-6630
Tel: (800) 543-8823
Fax: (856) 435-8182

www.biopharm.saint-gobain.com

FLS-5026C-0820-BPS

©2020 Saint-Gobain Life Sciences