

# Quick-Lock Connectors


## If it has to be quick!



**THE No. 1 System  
for Solvent Safety**

### Advantages

- ✓ **Drip protection**, thanks to integrated locking valve
- ✓ **Closed system**, for drip-free separation
- ✓ Very little air ingress
- ✓ **Minimal contamination**
- ✓ No outwashing of air bubbles from the system is necessary
- ✓ Intuitive connection
- ✓ User-friendly, ergonomic design enables a safe and reliable connection
- ✓ Chemically resistant materials
- ✓ Secure connection and reliable performance for many applications

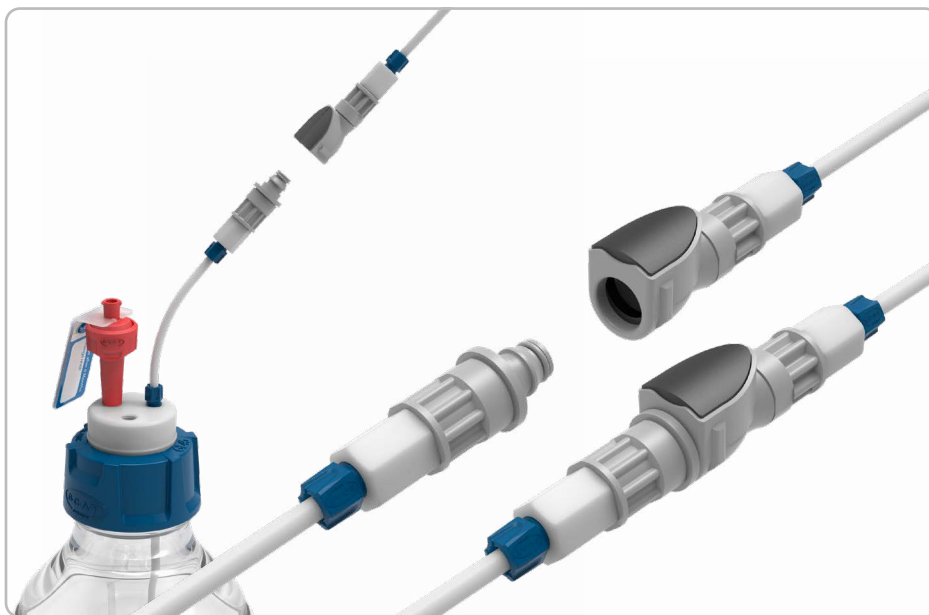

### Specifications

<b>Pressure</b>	Vacuum 15 psi @ 1.0 bar
<b>Temperature</b>	0°C to 49°C (32°F to 120°F)
<b>Materials</b>	<b>Main Components:</b> Glass-filled Polypropylene, PEEK, Polypropylene <b>Lock:</b> PEEK® <b>Springs:</b> Stainless Steel 316 <b>O-Rings:</b> EPDM
<b>Colour</b>	Grey
<b>Capillary Size:</b>	3.2 mm OD (1/8") 2.3 mm OD 1.6 mm OD (1/16")
<b>Discharge</b>	0.02 cc per disconnection @ 15 psi
<b>Uptake</b>	0.01 cc per connection

#### **A Quick-Lock Connector Set 1**

**360 189**



- For safe connecting of capillaries
- Suitable for capillaries with outer diameters of 3.2 / 2.3 / 1.6 mm
- PFA fittings in all 3 sizes are included in the scope of supply
- Contains Part Nos. 360 190 + 360 179

#### **B Quick-Lock Connector (m)**

**360 190**


- Spare parts for Set 360 189
- Incl. 3 PFA fittings (3.2 / 2.3 / 1.6 mm ID)



#### **C Quick-Lock Connector (f)**

**360 179**

- Spare parts for Set 360 189
- Incl. 3 PFA fittings (3.2 / 2.3 / 1.6 mm ID)

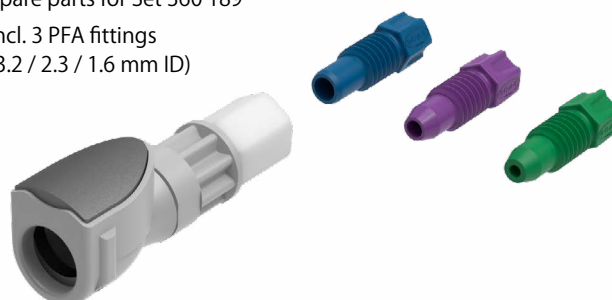



360 179

360 180

### D Quick-Lock Set II (for Safety Caps)

360 183


- Fast connecting to, and disconnecting from, HPLC units and solvent bottles, with a "click"!



- Suitable for capillaries with outer diameters 3.2 / 2.3 / 1.6 mm
- PFA fittings and sealing balls, in all 3 sizes, are contained within the scope of supply
- Consists of Part Nos. 360 180 + 360 179

### E 360 180


- Quick-lock connector (m) / Adapter for Safety Caps
- Incl. 3 PTFE sealing balls (3.2 / 2.3 / 1.6 mm ID)

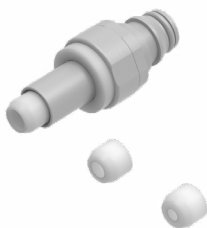

### C Quick-Lock Connector (f)

360 179

- Spare parts for Set 360 183
- Incl. 3 PFA fittings (3.2 / 2.3 / 1.6 mm ID)

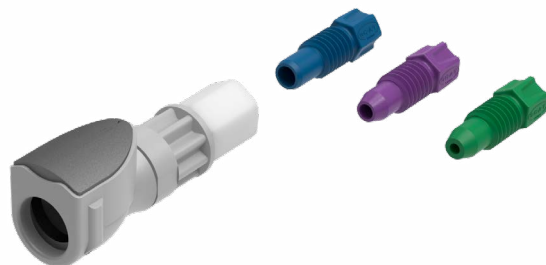

## Ordering Information

Fig.	Part No.	Description	PU
<b>A</b>	<b>360 189</b>	Quick-Lock Connector Set I for capillaries Contains 360 190 ( <b>B</b> ) and 360 179 ( <b>C</b> )	1 pc.
<b>B</b>	<b>360 190</b>	Quick-Lock Connector (m)	1 pc.
<b>C</b>	<b>360 179</b>	Quick-Lock Connector (f)	1 pc.
<b>D</b>	<b>360 183</b>	Quick-Lock Connector Set II for Safety Caps With adapter for screwing into the standard connectors of Safety Caps (UNF 1/4" 28G)	1 pc.
<b>E</b>	<b>360 180</b>	Quick-Lock Connector (m) for screwing into the standard connectors of Safety Caps (UNF 1/4" 28G)	1 pc.

### SCAT - Europe GmbH

Waldecker Straße 7  
64546 Mörfelden-Walldorf  
Germany

+49 (0) 61 05 / 30 55 86 - 0

+49 (0) 61 05 / 30 55 86 - 99

info@scat-europe.com

www.scat-europe.com

### SCAT - Americas, Inc.

1704 East Blvd., Suite 101  
Charlotte, NC 28203  
USA

+1.704.438.9477

+1.704.373.2603

info@scat-americas.com

www.scat-americas.com